



GLIDING-SKI-MODUL TYRAL



The two Gliding-Ski-Modules are also **easy to apply** in soft snow and allow the TYRAL stretcher to remain easy to slide and pull even in muddy snow.

The identical and symmetrical shape avoids confusing the attachment.



The Gliding-Ski-Modules can be attached in **three different** positions.

- For flat terrain, apply offset to the right and left (head and foot section).
- Attach to both sides of the foot for vertical downwards.
- Attach to both sides of the head for vertical upwards.

The two sliding ski modules can be easily attached to the stretcher **without tools** using a pipe clip.

Another advantage is that the two Glide-Ski-Modules can remain on the stretcher even when the wheel is in use - no constant conversion necessary. The wheel support is only removed for crossing snow fields.

A plastic tub was deliberately avoided in order to be able to offer a **sustainable, long-lasting** product.

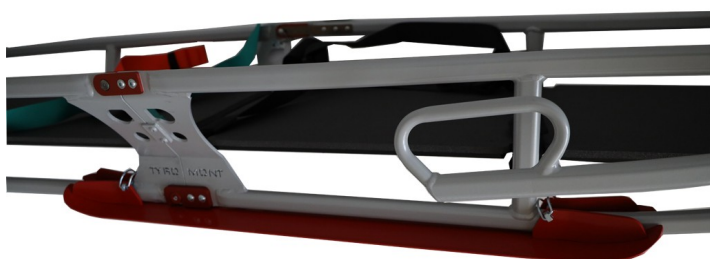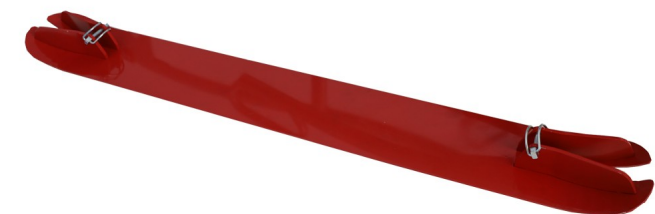

Art. No.	93262	
Dimensions	Length:	935 mm
	Width:	120 mm
Weight	Total per pair:	3,3 kg
Comments		



TYROMONT Alpin Technik GmbH
Bert-Köllensperger-Str. 6 · 6065 Thaur/Austria
T +43 (0) 5223 23 723-0 · office@tyromont.com
www.tyromont.com

Distributed by:
E-mail: sales@a-ss.net / www.a-ss.no / +47 3345 0910
Wirgenesvei 8, N-3157 Barkåker, NORWAY



Aviation and Survival Support
TRUST YOUR EQUIPMENT